

Children



IMPORTANT

- 1** In treating the heart area (4th intercostal, left of sternum) NEVER use a frequency other than 5 Hz.
- 2** coMra therapy is very good for children, especially as a form of preventative medication.
- 3** Any of the treatments prescribed for adults can be used on children as well, except that the time can be halved for infants under the age of 3 years, as infants and children respond much more quickly to laser treatment.
- 4** The few treatments given in this chapter are the ones that differ slightly from those given for adults. The rest are the same, the only difference being that you should always halve the exposure times given for adults when treating infants up to the age of 3 years.



Treatment times given in this user guide were designed for the home user with the Delta or Palm 980 nm laser to give a conservative starting point. Treatments with the professional 905 nm laser are more intense, and therefore these treatment times may be shortened, and also the space between treatments extended further, as determined by feeling and experience.




In treating the heart area NEVER use a frequency other than 5 Hz.

SECTION 17

PEDIATRICS 1

17 mins

Treatment time.

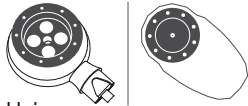

Using:
Delta Medical Terminal or
coMra Palm.

1 x daily

for 7 x days


• SOMATIC BIOSTIMULATION

TREATMENT POINT	⌚ MINS	⚡ HZ
1 Manubrium sterni	1 per point	5
2 4 th intercostal, left of sternum (only once per day)		
3 Tonsils		
4 Elbow fossae		
5 Right hypochondrium		
6 Epigastrium		
7 Left hypochondrium		
8 Scan along both sides of vertebral column	2 per side	1000
9 Kidneys	1 per point per frequency	5 then 50

i Somatic biostimulation is an excellent treatment for children older than 3 years who are weak and sickly. For very young infants and babies see MERIDIAN 2.



Treatments with the professional 905 nm laser are more intense, and therefore these treatment times may be shortened, and also the space between treatments extended further, as determined by feeling and experience.

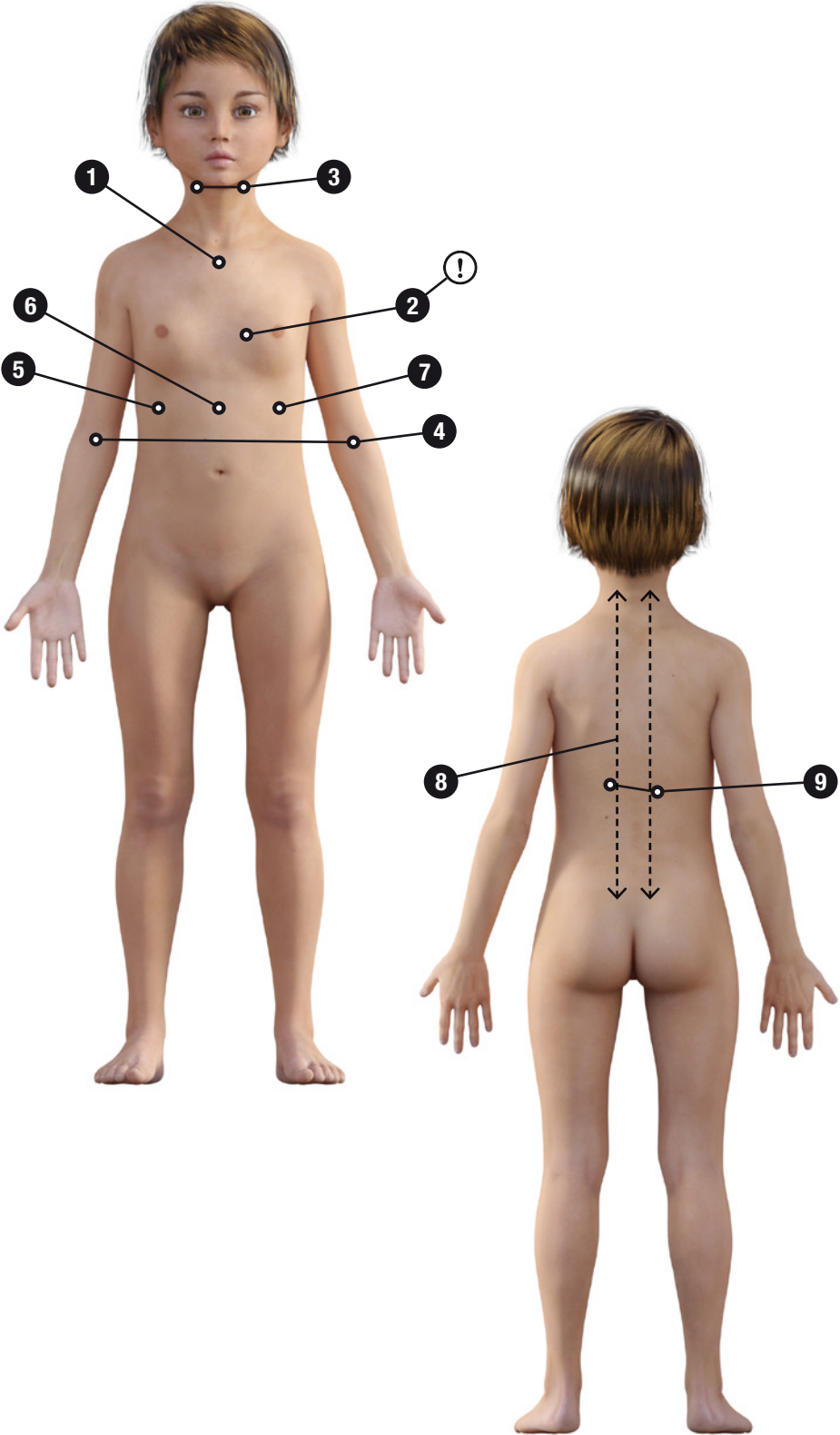


SECTION 17

PEDIATRICS 2

10 mins

Treatment time.

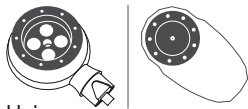

Using:
Delta Medical Terminal or
coMra Palm.

1 x daily for 10 x days


Thereafter every 2 months
if necessary.

• GASTROENTEROLOGY

TREATMENT POINT	⌚ MINS	⚡ HZ
1 Right iliac	1 per point	5
2 Right hypochondrium		
3 Epigastrium		
4 Left hypochondrium		
5 Left iliac		
6 Umbilicus		
7 Scan along both sides of vertebral column, mid thoracic region to lumbar	2 per side	1000



Treatments with the professional 905 nm laser are more intense, and therefore these treatment times may be shortened, and also the space between treatments extended further, as determined by feeling and experience.

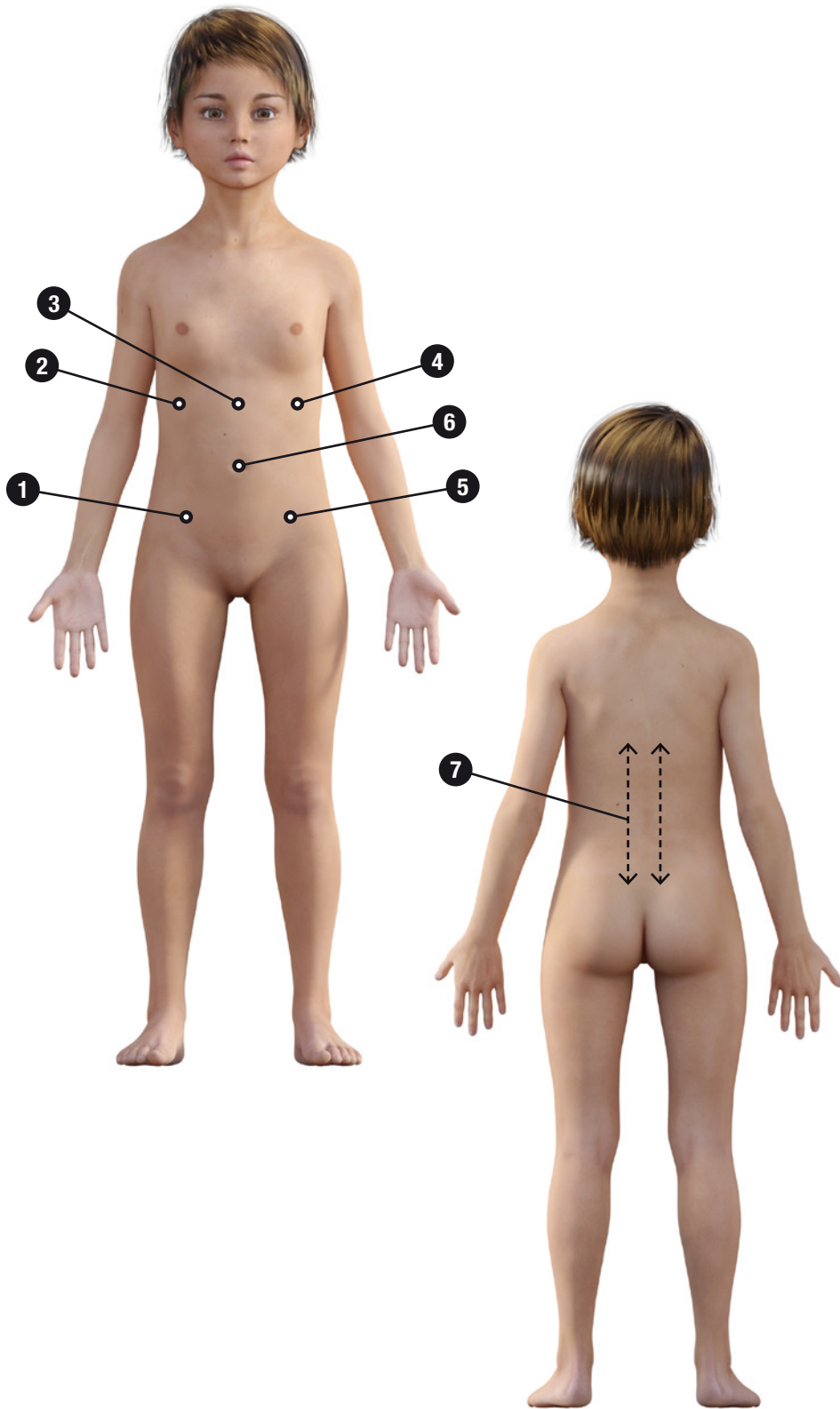


SECTION 17

PEDIATRICS 3

14 mins

Treatment time.

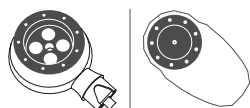

Using:
Delta Medical Terminal or
coMra Palm.

1 x daily for 10 x days

Thereafter every 2 months
if necessary.

• DISEASES OF THE KIDNEYS AND URINARY TRACTS

TREATMENT POINT	⌚ MINS	⚡ HZ
1 Bladder	2 per point	5
2 Bladder sphincter		
3 Kidneys		
4 Perineum		
5 Scan along both sides of vertebral column, lower back	2 per side	1000



Treatments with the professional 905 nm laser are more intense, and therefore these treatment times may be shortened, and also the space between treatments extended further, as determined by feeling and experience.

